

Comparison of DXA and Quantitative CT Bone Mineral Density Measurement using Contrast Enhanced Scans at the Proximal Femur: Implications for Opportunistic Osteoporosis Screening

Maciejewski A¹, Ziemelewick T¹, Binkley N¹, Brett A², Brown JK², Pickhardt PJ¹

¹University of Wisconsin School of Medicine and Public Health, Madison, WI, USA

²Mindways Software Inc, Austin, TX, USA

Introduction

We investigated IV contrast enhancement on femoral neck CTXA compared with DXA

For patients undergoing routine contrast-enhanced MDCT examinations, an opportunity exists for concurrent osteoporosis screening without additional radiation exposure or patient time using proximal femur quantitative CT (QCT) "CTXA". Previous studies demonstrated equivalence between unenhanced CT and DXA for femoral neck BMD evaluation; and have suggested a small BMD bias correction for enhanced CT studies. We investigated the effect of IV contrast enhancement on femoral neck CTXA T-score measurement compared with DXA.

Methods


Patients underwent contrast-enhanced CTXA and DXA

This cohort included 78 male and 277 female adults (mean age, 59.7±13.3 years; range, 21-90 years) who underwent standard contrast-enhanced CT assessment including the pelvis at 120kVp (GE Healthcare, Waukesha, WI) between November 2001 and December 2009. All subjects also had DXA BMD (GE Lunar, Madison, WI) assessment within 100 days of the CT study (mean 46±30 days). DXA T-scores were calculated using the NHANES III female reference data. Areal BMD in g/cm² of the femoral neck was measured on the CT series using QCT Pro Version 5.1 (Mindways Software, Austin, TX) with asynchronous phantom calibration. QCT T-scores were derived using manufacturer's female reference data.

Results

CTXA and DXA T-scores showed good correlation and no trend in their differences

Linear regression analysis showed good correlation between DXA and CTXA (R² = 0.824 for both BMD and T-scores) and the SD of the distribution of residuals (SEE) was 0.063 g/cm² or 0.45 T-score units. A Bland-Altman plot of T-scores indicates no trend in differences between the two measurements and a small bias with DXA T-score +0.18 units higher than CTXA.

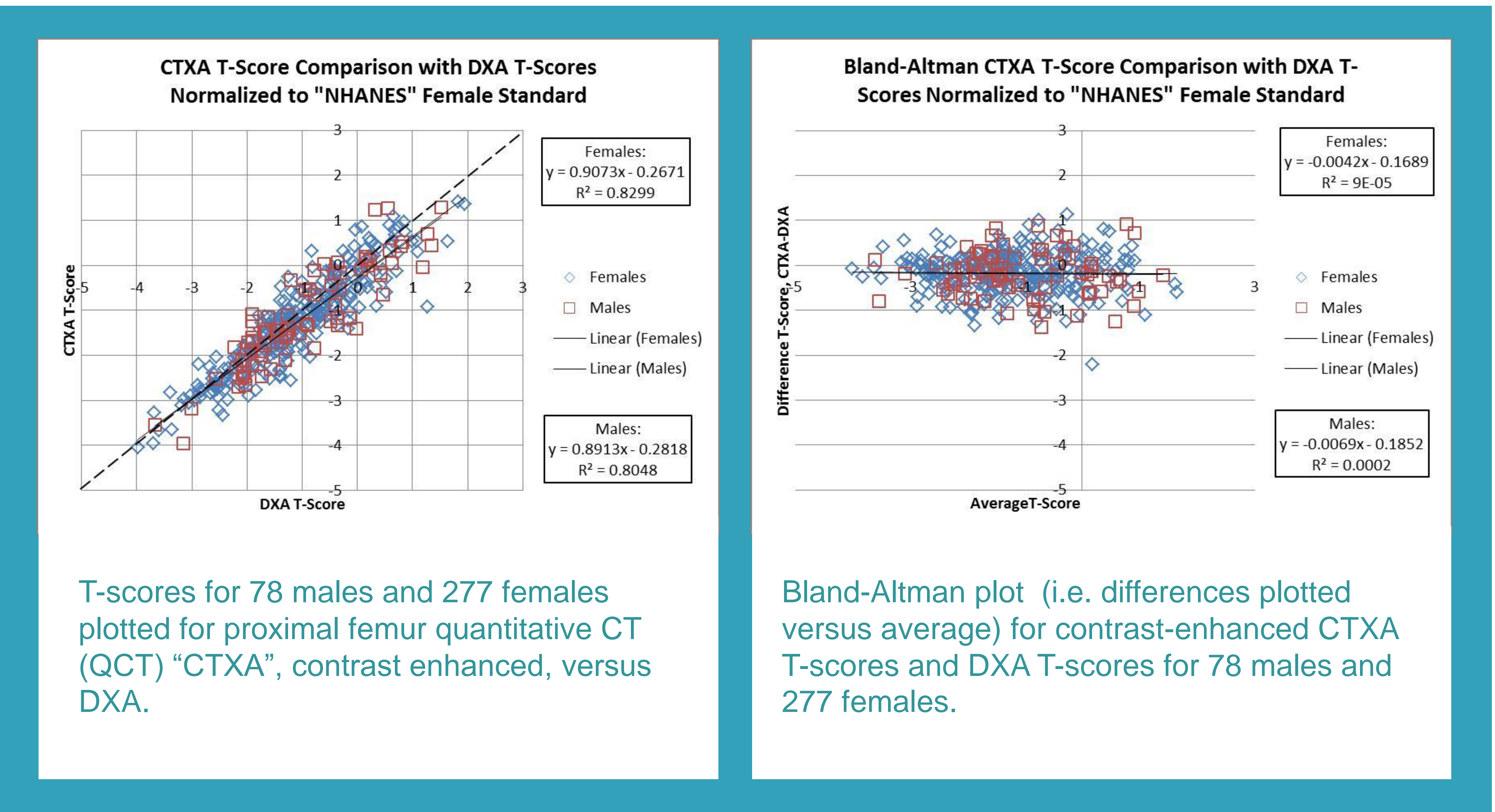

T-scores for 78 males and 277 females plotted for proximal femur quantitative CT (QCT) "CTXA", contrast enhanced, versus DXA.

Bland-Altman plot (i.e. differences plotted versus average) for contrast-enhanced CTXA T-scores and DXA T-scores for 78 males and 277 females.

Conclusion

Contrast-enhanced CTXA produces T-scores similar to DXA; now included in FRAX[®]

In this cohort, femoral neck T-scores obtained by CTXA and DXA are highly correlated. The SEE for BMD correlation compares favorably to that of a prior study comparing DXA and CTXA (0.053 g/cm²) and may be higher due to variance in bias introduced by contrast enhancement in this study. The small bias between DXA and CTXA T-scores presented here may be due to differences in CT scanner calibration due to differences in CT table scan height and scanner type. In conclusion, for opportunistic osteoporosis screening at routine post-contrast abdominopelvic CT scans, CTXA produces T-scores similar to DXA. Femoral neck QCT BMD measurement is now included in the WHO FRAX[®] tool and this method could greatly enhance osteoporosis screening since it can be applied regardless of the clinical indication for CT scanning.